



OPERATING INSTRUCTIONS MODEL R300 SERIES

SPEAKING


It is easy to control your speech with only one hand.

NO TOOLS REQUIRED.

Adjustable top head: Twist for fine-tuning of the sound, if needed.
WARNING - over adjusting of the head will result in diminished tone quality; be careful.

On/Off Button: To speak, push the button and release when you are finished.

Pitch Wheel: (on the left) Turn the wheel for easy adjustment to a higher or lower pitch. No need to further adjust after you have achieved your desired pitch level.

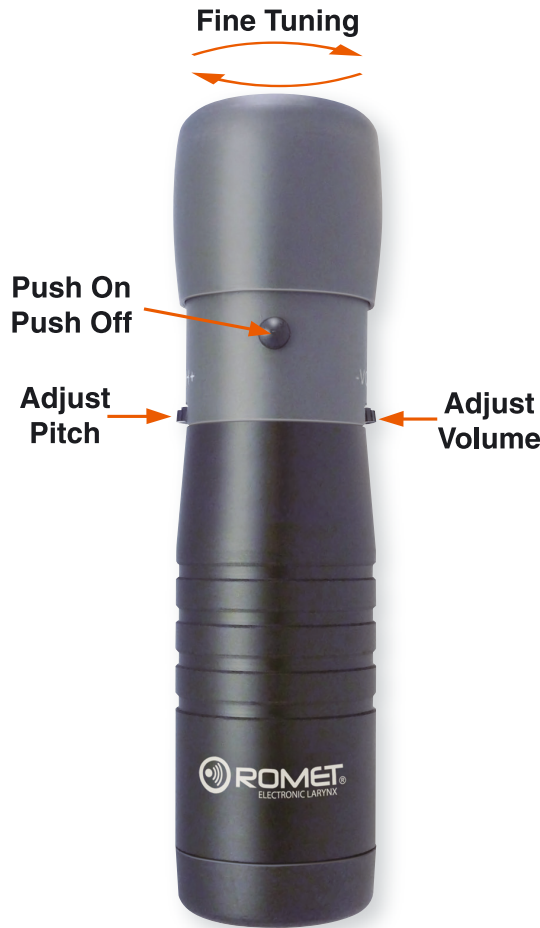

Volume Wheel: (on the right) Turn the wheel for easy adjustment while you are speaking to increase or decrease your speaking volume for emphasis.

Practice: Place the Romet® Electronic Larynx to your neck and practice finding your "sweet spot" to achieve your best vocal tone.

NOTE: This spot can be many places on either your neck or the side of your mouth; each person is different depending upon your skin texture and the amount of radiation you have received. Look into a mirror and repeat simple sentences to achieve vocal clarity. If possible, record your voice and play it back to hear how you sound. Practice regularly!

BATTERY INSERTION



Battery Insertion: Remove the battery end cap from the unit and insert the 9V Lithium battery into the barrel. The polarity of the battery is interchangeable; the + or - battery poles are reversible. Electrical contact will be made automatically when the barrel end cap is screwed onto the end of the barrel.

Slight pressure may be required the first time battery is inserted into unit. **DO NOT FORCE.**

BATTERY CHARGING



Charger: Included with your Model R300 unit you will find two 9V Lithium-Ion batteries. When you receive your unit the batteries are only approximately 80% filled to capacity. Charge your battery to full capacity prior to use to achieve maximum power. Plug your Lithium charger into a 110V circuit (220V


will require an adapter for your local area - NOT included in your kit). Insert one or two rechargeable batteries into your charger. Red LED light = Needs Charging. Green LED light = Fully Charged.

WARNING: Charge Li-Ion batteries in only a Li-Ion charger. DO NOT attempt to recharge them with any other charger. DO NOT put any other batteries in the charger, i.e. NiMH. **Normal 9V alkaline batteries will operate in your Romet® electro-larynx.**

ORAL ADAPTER INSTRUCTIONS



THE ORAL ADAPTER AND TUBES PROVIDED WITH YOUR ROMET ELECTRONIC LARYNX WILL ENABLE YOU TO BE ABLE TO BEGIN SPEAKING WHILE YOUR NECK IS HEALING POST SURGERY.

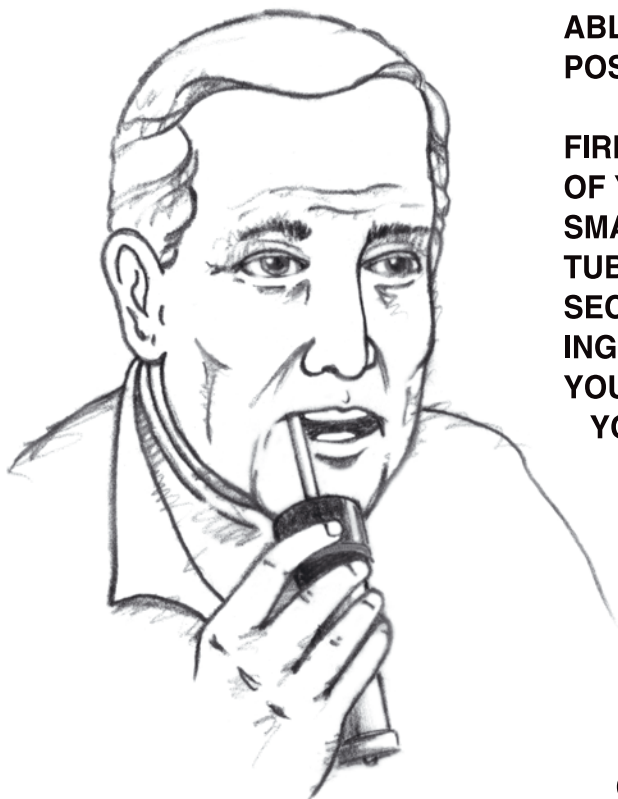

FIRMLY PLACE THE RUBBER CAP ADAPTER ON THE TOP OF YOUR ELECTRONIC LARYNX. INSERT A TUBE INTO THE SMALL HOLE ON THE TOP OF THE ADAPTER. PUSH THE TUBE FIRMLY INTO THE HOLE. THE TUBE MUST BE SECURE IN THE ADAPTER HOLE WHILE YOU ARE SPEAKING. PLACE THE END OF THE TUBE INTO THE CORNER OF YOUR MOUTH POINTING THE TIP TOWARDS THE TOP OF YOUR MOUTH.

ACTIVATE THE ON-OFF BUTTON. PRACTICE YOUR SPEECH WITH SIMPLE VOWELS OR SHORT SOUNDS. TRY A FEW SIMPLE SENTENCES.

IF POSSIBLE, SPEAK SLOWLY INTO A TAPE RECORDER. REPLAY YOUR SPEECH.

CLEAR TUBES ARE FITTED WITH TIPS TO PREVENT ANY FOOD PARTICLES FROM LODGING INTO THE TUBE WHILE BEING USED. YOU CAN EASILY REMOVE THIS TIP OR SHORTEN THE LENGTH OF THE TUBE WITH SCISSORS.

AFTER FINISHED SPEAKING, RINSE THE TUBES WITH COOL CLEAR WATER. BE SURE TO CLEAN THE TOP OF YOUR ELECTRONIC LARYNX TO REMOVE ANY SALIVA RESIDUE OR FOOD PARTICIPLES.



SANITAS GRUP Ltd.

1504 Sofia, Bulgaria, 25 "Tzar Osvoboditel" bul., ap. 1
office number: +359 2 943-38-38, office/fax: +359 2 944-16-50